


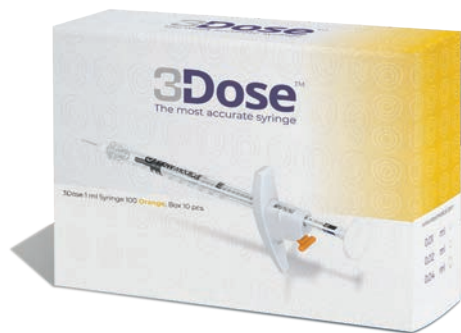

3DO100

3Dose 1 ml Syringe 100 Orange

Box with 10 syringes in sterile blister
(needles not included)

Settings:


0.01 ml - 0.02 ml - 0.04 ml

Used for reconstitution 4 U / 0.1 ml

Used for reconstitution 10 U Speywood / 0.1 ml

Used for reconstitution 20 U Speywood / 0.1 ml

Application: Wrinkle therapy

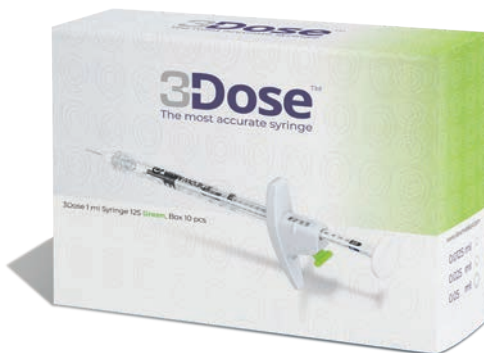

3DG125

3Dose 1 ml Syringe 125 Green

Box with 10 syringes in sterile blister
(needles not included)

Settings:

0.0125 ml - 0.025 ml - 0.05 ml

Used for reconstitution 5 U / 0.1 ml

Application: Wrinkle therapy, Mesotherapy


Dilution Table

Products	Reconstitution	Units / 0,1 ml	Plunger Setting	Number of units per click 1 click =
Neuromodulator Liquid (Galderma)	0,625 ml (ready to use, already reconstituted)	— 20 Speywood units / 0,1 ml	 green: 0,0125 — 2,5 black: 0,025 — 5 blue: 0,05 — 10	Speywood units Speywood units Speywood units
Neuromodulator Powder (Galderma)	0,625 ml (to be reconstituted)	— 20 Speywood units / 0,1 ml	 green: 0,0125 — 2,5 black: 0,025 — 5 blue: 0,05 — 10	Speywood units Speywood units Speywood units
	1,25 ml (to be reconstituted)	— 10 Speywood units / 0,1 ml	 green: 0,0125 — 1,25 black: 0,025 — 2,5 blue: 0,05 — 5	Speywood units Speywood units Speywood units
Neuromodulator Powder	50 U / 1,25 ml (to be reconstituted)	— 4 units / 0,1 ml	Green 3Dose  green: 0,0125 — 0,5 black: 0,025 — 1 blue: 0,05 — 2	units units units
Neuromodulator Powder	50 U / 1,0 ml (to be reconstituted)	— 5 units / 0,1 ml	Orange 3Dose  green: 0,01 — 0,5 black: 0,02 — 1 blue: 0,04 — 2	units units units

EN

PREDICTABLE[®] perfection

you can treat
40% more patients
with the
same amount of toxin
in a **shorter time**

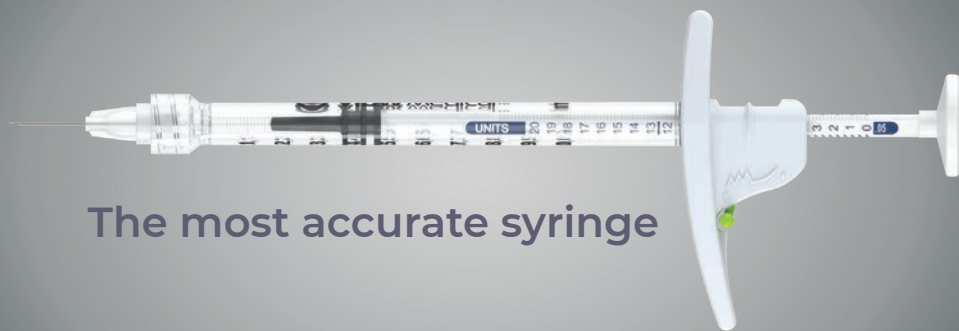

The most accurate syringe

"Once you start using the 3Dose syringes, you don't want to use anything else anymore"

Dr. Dalvi Humzah, Plastic Surgeon, DH Aesthetic Training

"We know that we do exactly and reproducibly, always get 1 unit of toxin per click!"

Dr. Sabine Zenker, Dr. Zenker Dermatology

"The 3Dose™ botulinum toxin syringe is one way to avoid errors"

Anna Baker, Award-Winning Aesthetic Nurse

"The syringe is easy to use and to read number of units for different dilutions"

Dr. Souphiyeh Samizadeh & Dr. Koenraad de Boulle

Fully focus
on the patient



Product Values



Saving Money

You can treat 40% more patients with the same amount of toxin in a shorter time

Working with the 3Dose™ syringe is more efficient, faster and the syringe has a low deadspace. By injecting just the right dose of toxin, overdose is a thing of the past. You can treat 40% more patients with the same amount of toxin in a shorter time. This accurate, optimized and faster working method allows you to achieve a significantly better return in your practice. Saving money, a good choice.



Accurate

You are assured of a predictable, accurate injection with every click

From now on you can fully focus on the patient and the treatment. You will find that the accuracy of the 3Dose™ syringe greatly improves the ease and efficiency of injections. You no longer need to calculate the units. During treatment, you can fully focus on where you want to inject, rather than how much to inject. You not only experience this convenience as very pleasant, but it also gives you a comfortable feeling.



Consistent

Your patient will be happier with consistent results

Because you inject smaller volumes more accurately, the volumes are exactly the same for each treatment (left/right, now/future). Because you work more precisely, the end result also looks better and more beautiful. This way you keep more satisfied patients because they notice that they get the best treatment and results every time.

User Benefits

Use less toxin & Exact Dosing

Low dead space save up to 20%. Predictable and reproducible, no over- and under-doses.

Focus

Tactile and audible feedback. Focus only on patient and case.

Calculation ml-unit

Is included in the syringe and setting.

Locking handle for

- Free moving-mode: plunger move free (pull up);
- Click-mode: operational, one click per dose.

Ergonomic handle

- Handle;
- Thumb-rest.

Display

- Remaining units;
- Stand upright.

